



MDF-DU900V-PE

VIP

-86°C Upright Freezer

845 L



PHCbi's MDF-DU900V **VIP** -86°C Upright Freezer provides maximum sample storage capacity within an optimum footprint combined with industry leading cooling performance and reliability for an ideal long-term preservation solution.

Maximum Sample Storage

The use of space saving patented VIP PLUS panels within the freezer cabinet allows a capacity of up to 672 2" boxes inside a footprint of just 1m² for maximum storage. The freezer's conventional depth allows for easy installation.

Reliable PHCbi Technologies

Cool Safe[®] compressors that are specifically designed for ultra-low temperature applications are employed in the proven PHCbi cascade refrigeration system ensuring the highest levels of performance and reliability.


Enhanced Use & Intelligent Security

The new EZlatch makes access to stored samples even easier. A colour LCD touch panel allows full user control, even with gloved hands, while the USB port makes transferring logged data to a PC simple and convenient.



High Capacity Sample Storage

Ideal for high volume and long term sample storage in biobanks and biorepositories.

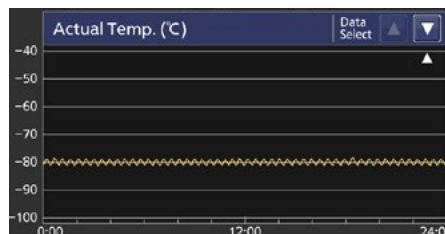

Dependable Preservation

Quality of design and manufacture ensures trusted and reliable storage that maintains the integrity of precious samples.



Easy Data Monitoring

Important information such as freezer temperature, door opening times and alarm history is logged for monitoring in GLP applications.

VIP -86°C Upright Freezer



Application Specific Design

Refrigeration components are carefully selected and matched for optimum and efficient operation under demanding laboratory environments. This improves temperature recovery after door openings and sample placement, while ensuring lower internal temperatures for increased durability.


Highest Box Capacity

The use of highly efficient VIP PLUS insulation optimizes interior volume in the smallest footprint possible. This ensures the most efficient use of space for up to 40%* more 2" box storage capacity within the same footprint.

*MDF-DU900V-PE compared with PHCbi model with an equivalent footprint.

Microprocessor Control with Touch Screen Display

Managed and monitored by an integrated microprocessor controller with a comprehensive alarm system and self- diagnostic functions. Freezer temperature and important status information are clearly displayed and logged, with the possibility to export the data to USB.

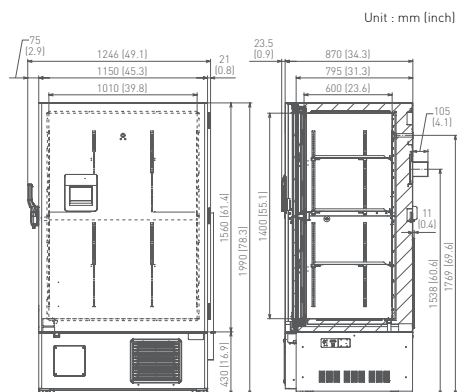

The MDF-DU900V-PE is certified as a Class IIa Medical Device (93/42/EEC and 2007/47/EC) for medical purposes of storing cells, tissues, organs and embryos.



PHC Europe B.V.
Nijverheidsweg 120 | 4879 AZ Etten-Leur | Netherlands
T: +31 (0) 76 543 3839 | F: +31 (0) 76 541 3732
www.phchd.com/eu/biomedical

Model Number		MDF-DU900V-PE
External Dimensions (W x D x H) ¹⁾	mm	1150 x 870 x 1990
Internal Dimensions (W x D x H)	mm	1010 x 600 x 1400
Volume	litres	845
Net Weight	kg	372
Capacity	2" boxes	672
Performance		
Cooling performance ²⁾	°C	-86
Temperature setting range	°C	-50 to -90
Temperature control range ²⁾	°C	-50 to -86
Control		
Controller		Microprocessor, non-volatile memory
Display		LCD Touch Screen
Temperature sensor		Pt-1000
Refrigeration		
Refrigeration system		Cascade
Compressors	W	2 x 1100
Refrigerant		HFC
Insulation material		PUF / VIP Plus
Insulation thickness	mm	70
Construction		
Exterior Material		Painted Steel
Interior Material		Painted Steel
Outer Door	qty	1
Outer Door Lock		Y
Inner Doors	qty	2 pieces (insulated)
Shelves	qty	3
Max. load - per shelf	kg	50
Max. load - total	kg	150
Vacuum release port		Y
Access Port	qty	2
Access Port Position		Back, Bottom
Access Port Diameter	Ø mm	17
Casters	qty	4 (2 leveling feet)
Alarms (R = Remote Alarm, V = Visual Alarm, B = Buzzer Alarm)		
Power Failure		V-B-R
High Temperature		V-B-R
Low Temperature		V-B-R
Filter		V-B
Door open		V-B
Electrical and Noise Level		
Power Supply	V	230
Frequency	Hz	50
Noise Level ³⁾	dB [A]	52
Options		
Small Inner Door Kit	set of 2	MDF-9ID-PW (max 2) ⁴⁾
Liquid CO ₂ back-up		MDF-UB6-PW
Temperature recorders		
- Circular type		MTR-G85C-PE ⁵⁾
- Chart paper		RP-G85-PW
- Ink pen		PG-R-PW
- Continuous strip type		MTR-85H-PW
- Chart paper		RP-85-PW
- Ink pen		DF-38FP-PW
- Recorder housing		MDF-S3085-PE
Optional communication systems		
Ethernet interface (LAN)		MTR-L03-PW

¹⁾ Exterior dimensions of main cabinet only, excluding handle and other external projections - See dimensions drawings on website for full details.

²⁾ Air temperature measured at freezer centre, ambient temperature +30°C, no load.

³⁾ Nominal value - Background noise 20dB(A).

⁴⁾ Installation of small inner door kit may effect usable storage capacity.

⁵⁾ Requires sensor cover MTR-DU700SF-PW.